MicroScan

MicroScan 0.5% N,N-Dimethylalphanaphthylamine Reagent Procedural Manual

REF B1010-45A, B1015-45

TABLE OF CONTENTS

English	2
Français (FR)	6
Deutsch (DE)	10
Italiano (IT)	14
Español (ES)	18
Português Portugal (PT-PT)	22
Norsk (NO)	26
Ελληνικά (EL)	30
中文 (ZH-CN)	34
Lietuviškai (LT)	37
Magyar (HU)	41
Polski (PL)	45
Čeština (CZ)	49
Slovenčina (SK)	53
한국어 (KO)	57
Türkçe (TR)	61
Русский (RU)	65
Hrvatski (HR)	70
Български (BG)	74
中文 (ZH-TW)	78
Română (RO)	81
Србија (SR)	85
Português Brasil (PT-BR)	89
Nederlands (NL)	93
Tiếng Việt (VI)	97
Македонски (MK)	101

MicroScan

MicroScan 0.5% N,N-Dimethylalphanaphthylamine Reagent Procedural Manual

REF B1010-45A, B1015-45

INTENDED USE

MicroScan 0.5% N,N-Dimethylalphanaphthylamine Reagent is for use with MicroScan panels.

INTENDED USER

This product is intended for laboratory professional use.

CLINICAL RELEVANCE

Not applicable, product is an IVD accessory.

SUMMARY AND PRINCIPLE

MicroScan 0.5% N,N-Dimethylalphanaphthylamine is an *in vitro* diagnostic accessory which aids in the detection of nitrate reduction for qualitative bacterial identification. Refer to MicroScan Gram Neg ID Panel Type 2, Gram Pos ID Panel Type 3, HNID, and Rapid Anaerobe ID panel instructions for methods of incubation and results determination which may include manual, semi-automated, or full automation on MicroScan instrumentation.

Nitrate: The ability of an organism to reduce nitrate to nitrite is detected by the addition of 0.8% Sulfanilic Acid and 0.5% N,N-Dimethyl-alpha-naphthylamine, which produce a red color in the presence of nitrite.

Nitrite: The ability of an organism to reduce free nitrite to free nitrogen gas is detected by the addition of 0.8% Sulfanilic Acid and 0.5% N,N-Dimethyl-alpha-naphthylamine. In the presence of nitrite a red color develops which is a negative test, if no color develops, the nitrite has been reduced and the test is a positive test.

WARNING AND PRECAUTIONS

1. For *in vitro* diagnostic use only.
2. Observe aseptic techniques and established precautions against microbiological hazards throughout all procedures, with special awareness that the inoculated panels contain potentially pathogenic organisms.
3. Contains: N,N-Dimethylalphanaphthylamine and Glacial Acetic Acid.
4. Discard as appropriate for chemical waste.

STORAGE

MicroScan 0.5% N,N-Dimethylalphanaphthylamine Reagent must be stored at 15 - 30°C away from light (e.g. in dark amber glass containers). Do not use beyond the expiration date.

SHELF LIFE/EXPIRATION DATE

Expiration Date: The reagent is stable for one year after the date of manufacture.

Shelf Life: The reagent is stable for one year not to exceed expiration date on bottle when stored properly.

On-board Instrument Shelf Life: The reagent is stable on the instrument for two months not to exceed expiration date on bottle.

GHS HAZARD CLASSIFICATION

MicroScan 0.5%

N,N-Dimethylalphanaphthylamine

DANGER

H314

Causes severe skin burns and eye damage.

P280

Wear protective gloves, protective clothing and eye/face protection.

P301+P330+P331

IF SWALLOWED: rinse mouth. Do NOT induce vomiting.

P303+P361+P353

IF ON SKIN (or hair): Rinse skin with water.

P305+P351+P338

IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.
Immediately call a POISON CENTER or doctor/physician.
Acetic Acid 25 - 35%

P310

Safety Data Sheet is available at beckmancoulter.com/techdocs

EVIDENCE OF DETERIORATION

MicroScan 0.5% N,N-Dimethylalphanaphthylamine should be a colorless, clear, pungent liquid. Do not use beyond the expiration date. Reagent should not be used if packaging is punctured, unsealed or torn. Contact your Beckman Coulter Representative or Distributor for further assistance.

SPECIMEN COLLECTION AND PREPARATION

Refer to the appropriate MicroScan panel Procedural Manual.

MATERIALS PROVIDED

MicroScan 0.5% N,N-Dimethylalphanaphthylamine, catalog number B1010-45A, contains 30 mL

MicroScan 0.5% N,N-Dimethylalphanaphthylamine, catalog number B1015-45, contains 250 mL

COMPOSITION

N,N-Dimethylalphanaphthylamine 0.30 mL/mL

Glacial Acetic Acid 0.00516 mL/mL

This product does not contain materials from animal content.

MATERIALS REQUIRED, BUT NOT PROVIDED

Refer to the appropriate MicroScan panel Procedural Manual.

PROCEDURE

Following panel incubation, add 1 drop of 0.8% Sulfanilic Acid then 1 drop of 0.5% N,N-Dimethylalphanaphthylamine. Refer to the appropriate MicroScan Panel Procedural Manual for details on reaction and read time. The test can be performed either manually or using automation with WalkAway Systems (refer to the appropriate instrument Operators Guide for details).

RESULTS

Refer to the appropriate MicroScan panel Procedural Manual.

LIMITATION OF THE PROCEDURE

Refer to the appropriate MicroScan panel Procedural Manual.

QUALITY CONTROL

Refer to the appropriate MicroScan panel Procedural Manual.

PERFORMANCE CHARACTERISTICS

Refer to the appropriate MicroScan panel Procedural Manual.

REPRODUCIBILITY

Refer to the appropriate MicroScan panel Procedural Manual.

WARRANTY STATEMENT

The system is covered by and subject to the provisions of the warranty included in your contractual agreement for the system or its reagents. The customer is responsible for routine preventive maintenance procedures. Repairs arising from the failure to perform these maintenance procedures at the indicated time intervals are made at the discretion of Beckman Coulter, and at the customer's expense.

NOTICE TO USER

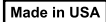
Any serious incident that has occurred in relation to MicroScan panels shall be reported to Beckman Coulter and the competent authority of the Member State in which the user and/or the patient is established.

REVISION HISTORY

REVISION	DATE	DESCRIPTION
AA	November 2021	Initial Release
AB	April 2022	Added: CE Mark

Symbols Key

The following symbols may be used in the labeling of Beckman Coulter Microbiology products.

Symbols Key			
Symbol	Symbol Title	Description of Standard	Standard
	Do not re-use	Indicates a medical device that is intended for one use, or for use on a single patient during a single procedure.	ISO 15223-1, clause 5.4.2
	Use-by-date	Indicates the date after which the medical device is not to be used.	ISO 15223-1, clause 5.1.4
	Batch Code	Indicates the manufacturer's batch code so that the batch or lot can be identified.	ISO 15223-1, clause 5.1.5
	Catalogue Number	Indicates the manufacturer's catalogue number so that the medical device can be identified.	ISO 15223-1, clause 5.1.6
	Manufacturer	Indicates the medical device manufacturer, as defined in EU Directives 90/385/EEC, 93/42/EEC and 98/79/EC.	ISO 15223-1, clause 5.1.1
	Date of Manufacture	Indicates the date when the medical device was manufactured.	ISO 15223-1, clause 5.1.3
	Authorized Representative in the European Community	Indicates the authorized representative in the European Community.	ISO 15223-1, clause 5.1.2
	Contains sufficient for <n> tests	Indicates the total number of IVD tests that can be performed with the IVD. (Typically included on reagent kits).	ISO 15223-1, clause 5.5.5
	<i>In Vitro</i> Diagnostic Medical Device	Indicates a medical device that is intended to be used as an <i>in vitro</i> diagnostic medical device.	ISO 15223-1, clause 5.5.1
	Temperature Limitation	Indicates the temperature limits to which the medical device can be safely exposed.	ISO 15223-1, clause 5.3.7
	Consult Instructions for Use	Indicates the need for the user to consult the instructions for use.	ISO 15223-1, clause 5.4.3
	CE Mark	European Conformity Regulatory Mark	n/a
	Contents	n/a	n/a
	Contents (package)	Antimicrobial Agent (abbreviation)	n/a
	Contents (package)	Identification Substrate (abbreviation)	n/a
	Safety Data Sheet	Indicates a safety data sheet.	n/a
	Made in USA	n/a	n/a
	For Export Use Only	n/a	n/a
	Reconstitution Volume	n/a	n/a
	Fragile, handle with care	Indicates a medical device that can be broken or damaged if not handled carefully.	ISO 15223-1, clause 5.3.1
	This way up	Indicates the correct upright position	ISO 7000: 0623
	Info For USA Only	Indicates information that is required only for products intended to be sold in the US.	n/a
	Prescription Use Only	Recognized by the US FDA as an alternate to "Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner"	21CFR 801.109(b)(1)

Symbols Key			
Symbol	Symbol Title	Description of Standard	Standard
	Biological Risks	Indicates that there are potential biological risks associated with the medical device. Note: This symbol is not to be confused with the "biohazard" sign intended to be used in the workplace.	ISO 15223-1, clause 5.4.1
	Danger	n/a	n/a
ISO 15223-1: Medical Devices - Symbols to be Used with Medical Device Labels, Labeling and Information to be Supplied. Part 1: General Requirements.			

Beckman Coulter, the stylized logo, and the Beckman Coulter product and service marks mentioned herein are trademarks or registered trademarks of Beckman Coulter, Inc. in the United States and other countries.

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Beckman Coulter Ireland Inc.,
Beckman Coulter, Inc., 250 S. Kraemer Blvd., Brea, CA 92821 U.S.A.